

Simplifying DXEK/DSAEK

Innovators



Prof. Donald Tan
Singapore



Tan DSAEK Forceps: AE-4226

"These forceps enable quite smooth insertion of the donor lamella into the anterior chamber. The curve and tip of the forceps is well designed, and suited for pull-through technique of DSAEK."

Akira Kobayashi, MD
Ishikawa, Japan

"It is just perfect!"

Prof. Paolo Rama
Milan, Italy



Thomas John, MD
Chicago, IL

John DXEK/DSAEK Double-Ended Marker: AE-2712

- Used to mark the patient cornea at either 8mm or 9mm based on the surgeon's preference



John Dexatome DXEK/DSAEK Spatula: AE-2872

- Single instrument allows for both Descemetorhexis and Descemet's membrane detachment and subsequent removal
- Unique multi-curve design facilitate complete, 360° scoring of Descemet's membrane without exiting the anterior chamber

DVD
VIDEO
AVAILABLE



Sinskey No Hole IOL Manipulator: AE-2209

- Modified Sinskey Hook designed to score the Descemet's membrane




John DXEK/DSAEK Descemet's Stripper: AE-2874

- The T-bar shape helps remove the Descemet's membrane as a single disk and prevents or limits uneven segmental tearing of Descemet's membrane during its detachment from the host inner corneal dome
- Curvilinear design is effective in reaching Descemet's membrane over all inner curvatures of the recipient cornea

DVD
VIDEO
AVAILABLE



John DXEK Inserting Forceps: AE-4227

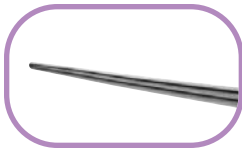

DVD
VIDEO
AVAILABLE

- Highly polished jaws are ideal for grasping and removing the damaged endothelium

John DXEK/DSAEK Stromal Scrubber: AE-2878



DVD
VIDEO
AVAILABLE

- Effectively roughens the inner corneal stroma within the Descemetorhexis circle
- Curvilinear design allows for easy access of the inner-corneal dome


Tan DSAEK Forceps: AE-4226



DVD
VIDEO
AVAILABLE

- Safely and effectively inserts the donor button with the pull-through technique using a lens glide
- Atraumatic ridges secure grip of stromal layer in the presence of visco

John DXEK/DSAEK Fixation Hook: AE-2182

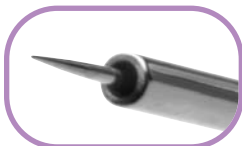

DVD
VIDEO
AVAILABLE

- Holds the donor button in place during air injection and retracts to prevent accidental drag when exiting the chamber

John DXEK/DSAEK Glider: AE-2879



DVD
VIDEO
AVAILABLE

- Designed to gently smoothen out wrinkles at the donor-host interface when used along the outer corneal surface