


# John DMEK Set

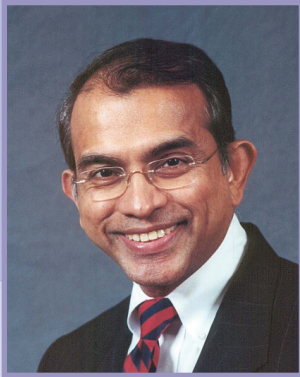

Dr. Thomas John  
Chicago, IL

***DMEK is a surgically challenging technique for the corneal surgeon. However, it is very gratifying to both the patient and the surgeon with regard to the level and quality of vision. DMEK has several distinctive advantages compared to DSAEK, including the restoration of the recipient corneal anatomy to near normal corneal structural configuration. Unlike DSAEK, which is an additive procedure, DMEK is to be considered as a substitution procedure that appears to be superior to DSAEK.***

## John DMEK/DSAEK Dexatome Spatula: AE-2872



To be used with the John DMEK Large Crystal Spatula for simultaneous unravelling of the donor DM within the recipient anterior chamber.

## John DMEK Twin Block: AE-1586



This Block consists of twin wells that is specially designed to hold the donor cornea and help the surgeon in the donor tissue preparation. It facilitates separation of the Descemet's membrane (DM) from the donor corneal button. The John DMEK block helps in separation, staining, and temporary placement of the donor DM and endothelial cells until transplanted to the recipient cornea.


## John DMEK Large Crystal Spatula: AE-2328



This spatula has a smooth, rounded, olive tip extremity which is specially designed for use on the donor Descemet's membrane within the recipient anterior chamber. It facilitates unrolling of the donor DM.

**Also available: John DMEK Small Crystal Spatula: AE-2329**

## John DMEK Cannula, 30g: AE-7003

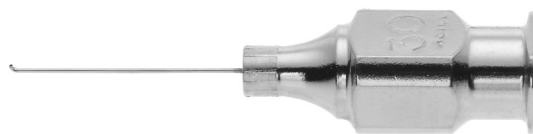

Facilitates irrigation of BSS and air injection into the recipient anterior chamber as needed during the transplantation of the DM to the recipient cornea.

*The more you know...*

**ASICO**<sup>®</sup>

Today's Precision... Tomorrow's Vision<sup>™</sup>

ASICO LLC

26 Plaza Drive • Westmont, IL 60559 • USA

Phone: 630-986-8032 • Fax: 630-986-0065

Customer Service: 1-800-628-2879

www.asico.com • E-mail: info@asico.com

© 2010 ASICO LLC



DVD  
VIDEO  
AVAILABLE