

ASICO[®] Toric Solutions

Problem


Solution

Feedback

Pre-Operative Marking

Pre-operative markings do not align to the horizontal axis.

Nuijts/Lane Pre-Op Toric Reference Marker: AE-2791TBL

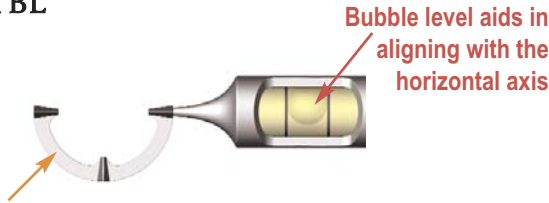

Non-marking surface is minimized to maximize visibility


“With the bubble level marker, I no longer need the slit lamp to mark the horizontal axis.”

Erik Mertens, MD, Antwerp, Belgium

Pre-operative markings do not last through surgery.

Pre-Op Toric Reference Marker (Single-Side): AE-2793S

Handle position allows pressure to be evenly distributed through out the marking points

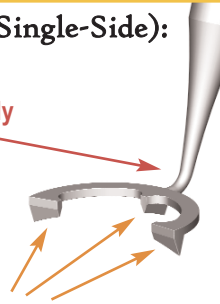

Marking surfaces contour to the curvature of the globe for accurate placement and effective marking

“This new design allows me to apply less pressure to make good marks in a manner that is non-threatening to my patients.”


Stephen S. Lane, MD, Minneapolis, MN

Intra-Operative Marking

Intra-Op marker does not align to the pre-operative markings.

Intra-Op Toric Axis Marker II: AE-2794

Cutouts at 0° and 180° allow you to identify and accurately align to the pre-operative reference marks

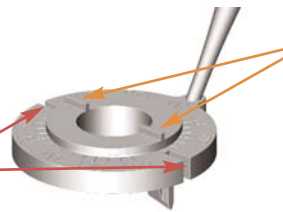

Dial references marking position and allows for easy rotation to the proper degree setting with the surgeon's fingers

“With this marker I can clearly see the pre-operative marks which allow me to accurately mark the alignment reference lines.”


David M. Lubeck, MD, Homewood, IL

Intra-Op marker obscures the pre-operative markings and is too large for patients with smaller eyes.

Modi Intra-Op Toric Marker II: AE-2848N

Lines indicate location of axis marking pattern

Wide opening for increased visibility and for centering corneal light reflex

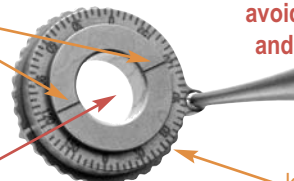

Small, 12.5mm diameter avoids obstruction and will not strain small eyes

Outer barrel knurled for easy axis rotation

“This marker allows the pre-operative markings to be seen easily for consistent alignment and works well on smaller corneas.”

Satish S. Modi, MD, Poughkeepsie, NY

The more you know...
ASICO[®]
Today's Precision... Tomorrow's Vision[™]

ASICO LLC
26 Plaza Drive • Westmont, IL 60559 • USA
Phone: 630-986-8032 • Fax: 630-986-0065
Customer Service: 1-800-628-2879
www.asico.com • E-mail: info@asico.com

© 2009 ASICO LLC

